


EXACTRAC 6.0 FOR ELEKTA™ TREATMENT SYSTEMS

Accurate, fast and intuitive IGRT for true frameless radiosurgery.

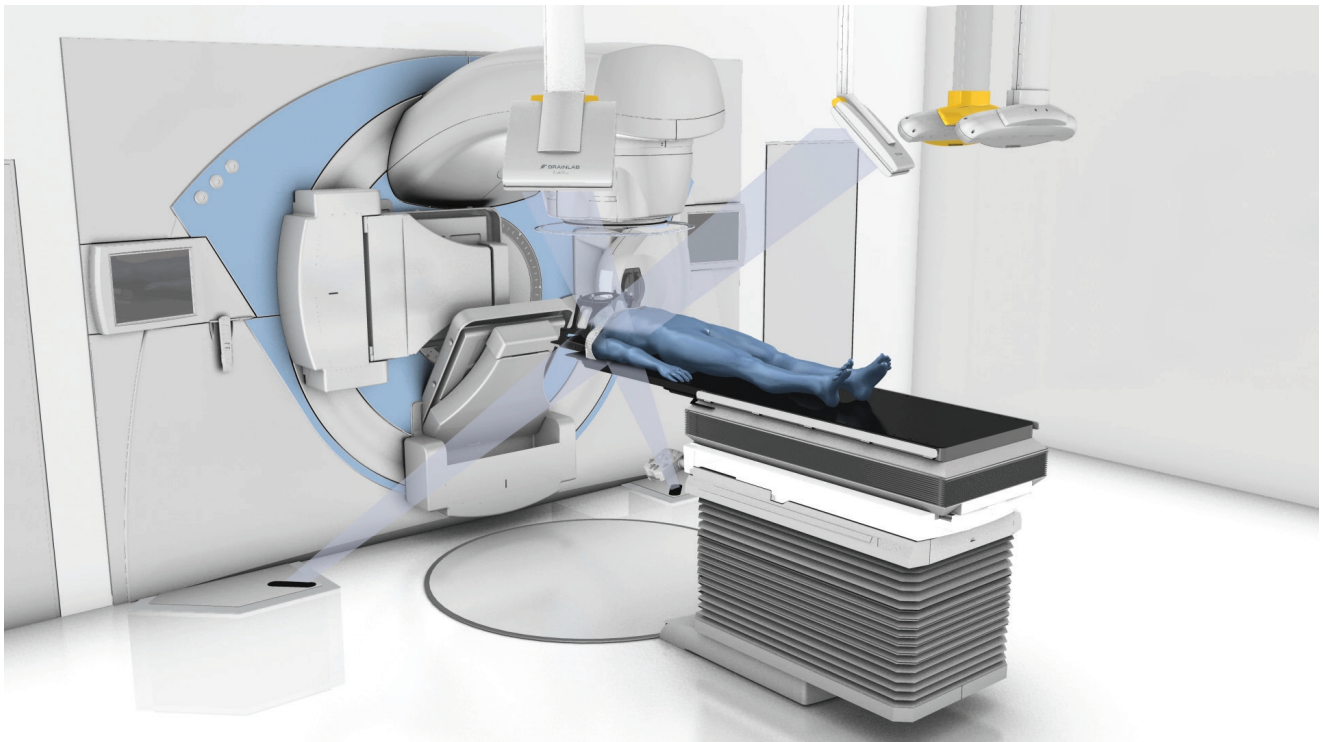

STREAMLINED COMPATIBILITY

- › The Exactrac® IGRT system provides a consistent and intuitive integration for Elekta linacs: Versa HD™, Axesse™, Infinity™, Precise Treatment System™ and Synergy®.

ENHANCED WORKFLOW

- › The room-based setup, combined with an integrated optical tracking system, allows for an efficient 2-minute patient setup.

FRAMELESS SRS TREATMENTS

- › Highly accurate delivery of SRS without conventional head ring based on non-invasive stereotactic mask system and robotic couch HexaPOD evo RT.

6DOF PATIENT ALIGNMENT WITH HEXAPOD EVO RT SYSTEM

- › Proprietary 6D fusion and robotic alignment for accurate 6DOF patient positioning.

HIGH ACCURACY

- › The calibration process of Exactrac, directly correlated to the radiation isocenter, contributes to a high intrinsic accuracy.

X-RAY MONITORING FOR INSTANT VERIFICATION DURING TREATMENT

- › Exactrac allows for instantaneous imaging of internal structures anytime during beam-on, even at various table and gantry angles.

NON COPLANAR POSITIONING VERIFICATION

- › With its capability to image internal structures at any table angle, Exactrac provides the confidence to maintain the precision of the initial setup throughout the treatment session.